


GSK343

## Chemical Properties

|                   |  |
|-------------------|--|
| CAS No. :         | 1346704-33-3   |
| Formula:          | C31H39N7O2   |
| Molecular Weight: | 541.69   |
| Storage:          | keep away from direct sunlight<br>Powder: -20°C for 3 years   In solvent: -80°C for 1 year |

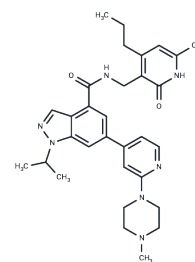

## Biological Description

|               |   |
|---------------|---|
| Description   | GSK343, a specific and effective EZH2 inhibitor (IC50=4 nM), exhibits 60 fold specificity activity against EZH1, and >1000 fold specificity activity against other histone methyltransferases.  |
| Targets(IC50) | Histone Methyltransferase, Autophagy  |
| In vivo       | Mice treated with GSK343 (5 mg/kg) displayed significantly inhibited tumor growth compared to the control group. The average volume and weight of tumors in the GSK343-treated group were notably decreased. This reduction in tumor growth was evident as early as 20 days post-implantation and continued throughout the study. Furthermore, in the GSK343-treated animals within the xenograft model, there was a substantial increase in the levels of E-cadherin mRNA, alongside a significant decrease in vimentin mRNA levels, compared to the control group[3].   |
| Kinase Assay  | In vitro biochemical assays against histone methyltransferases: Activity against EZH2 is assessed using 5 member PRC2 complex (Flag-EZH2, EED, SUZ12, AEBP2, RbAp48). The assay protocol may be summarized as follows: 10 mM stocks of compounds are prepared from solid in 100% DMSO. An 11 point serial dilution master plate is prepared in 384 well format (1:3 dilution, columns 6 and 18 were equal volume DMSO controls) and dispensed to assay ready plates using acoustic dispensing technology to create a 100 nL stamp of compound and DMSO controls. The assay additions consisted of equal volume additions of 10 nM EZH2 and the substrate solution (5 µg/mL HeLa nucleosomes and 0.25 µM [3H]-SAM) dispensed into assay plates using a multi-drop combi dispense. Reaction plates are incubated for 1 hr and quenched with an equal volume addition of 0.5 mg/mL PS-PEI Imaging Beads (RPNQ0098) containing 0.1 mM unlabeled SAM. The plates are sealed, dark adapted for 30 minutes, and a 5 minute endpoint luminescence image is acquired using a Viewlux imager. Plate statistics such as Z' and signal to background as well as dose response curves are analyzed using Activity BaseXE. The in vitro biochemical activity of EZH1 is assessed as part of a 5 member PRC2 complex using a 384 well SPA assay identical to EZH2. Buffer components, reagent dispensing, compound plate preparation, quench conditions and data analysis are identical for EZH1 and EZH2 with final assay concentrations of 20 nM EZH1, 5 µM/mL HeLa nucleosomes and 0.25 µM [3H]-SAM. Further data analysis, pIC50 pivots and visualizations are enabled by TIBCO Spotfire. Compounds are profiled at Reaction Biology Corp. (Malvern, PA) to assess inhibition in their panel of histone methyltransferase assays. Methyltransferase activity is assessed using HotSpot |

technology, a miniaturized radioisotope-based filter binding assay. Inhibitors are dissolved in dimethyl sulfoxide (DMSO) and tested at concentrations up to 100  $\mu$ M with a final DMSO concentration of 2%. Buffer containing the methyltransferase at the listed concentration and its preferred substrate as shown in the accompanying table is preincubated with compound for 10 min. Reactions are initiated by the addition of 1  $\mu$ M S-adenosyl-L-[methyl- $^3$ H]methionine (SAM), allowed to incubate for 60 min at 30C followed by transfer to P81 filter-paper and PBS wash before detection.

## Cell Research

To account for varying doubling rates among cancer cell lines, the optimal cell seeding is determined empirically for all cell lines by examining their growth in a 384-well plate over 6 days with a wide range of seeding densities. Cells are then plated at the optimal seeding density and allowed to adhere overnight. Cells are treated in duplicate with a 20-point 2-fold dilution series of compound or 0.147% DMSO (vehicle control) and incubated for 6 days at 37C in 5% CO<sub>2</sub>. Cells are then lysed with 25  $\mu$ l CellTiter-Glo per well and chemiluminescence is quantified with a TECAN Safire2 microplate reader. In addition, an untreated plate of cells is harvested at the time of compound addition (T<sub>0</sub>) to quantify the starting number of cells. CTG values after 6 days of treatment were expressed as a percent of the T<sub>0</sub> value and plotted against compound concentration. Data are fit with a 4-parameter equation to generate a concentration response curve and the concentration of compound required to inhibit 50% of growth (GI<sub>50</sub>) is determined (Only for Reference)

## Solubility Information

## Solubility

DMSO: 22.27 mg/mL (41.11 mM), Heating is recommended.  
 10% DMSO+40% PEG300+5% Tween 80+45% Saline: 1 mg/mL (1.85 mM), Sonication is recommended.  
 (< 1 mg/ml refers to the product slightly soluble or insoluble)

## Preparing Stock Solutions

|       | 1mg       | 5mg       | 10mg       |
|-------|-----------|-----------|------------|
| 1 mM  | 1.8461 mL | 9.2304 mL | 18.4607 mL |
| 5 mM  | 0.3692 mL | 1.8461 mL | 3.6921 mL  |
| 10 mM | 0.1846 mL | 0.923 mL  | 1.8461 mL  |
| 50 mM | 0.0369 mL | 0.1846 mL | 0.3692 mL  |

Please select the appropriate solvent to prepare the stock solution, according to the solubility of the product in different solvents. Please use it as soon as possible.

### Reference

Verma SK, et al. ACS Med Chem Lett. 2012, 3, 1091-1096.

Zhao Deng, Chenbin Cui, Yanan Wang, Jiangjin Ni, et al. FSGHF3 and peptides, prepared from fish skin gelatin, exert a protective effect on DSS-induced colitis via the Nrf2 pathway. Food & Function. 2020

Zhao Deng, Chenbin Cui, Yanan Wang, Jiangjin Ni, et al. FSGHF3 and peptides, prepared from fish skin gelatin, exert a protective effect on DSS-induced colitis via the Nrf2 pathway. Food & Function. 2020

Amatangelo MD, et al. Cell Cycle. 2013, 12(13), 2113-2119.

Wang Z, Guo Z, Wang X, et al. Inhibition of EZH2 Ameliorates Sepsis Acute Lung Injury (SALI) and Non-Small-Cell Lung Cancer (NSCLC) Proliferation through the PD-L1 Pathway. Cells. 2022, 11(24): 3958.

Ding M, et al. The polycomb group protein enhancer of zeste 2 is a novel therapeutic target for cervical cancer. Clin Exp Pharmacol Physiol. 2015 May;42(5):458-64.

Xing Y, Tan C, Liu Z, et al. Resveratrol as a BCL6 natural inhibitor suppresses germinal center derived Non-Hodgkin lymphoma cells growth. Journal of Natural Medicines. 2025: 1-13.

Quan C, Chen Y, Wang X, et al. Loss of histone lysine methyltransferase EZH2 confers resistance to tyrosine kinase inhibitors in non-small cell lung cancer. cancer letters. 2020

**Inhibitor · Natural Compounds · Compound Libraries · Recombinant Proteins**

This product is for Research Use Only · Not for Human or Veterinary or Therapeutic Use

Tel: 781-999-4286    E\_mail: info@targetmol.com    Address: 34 Washington Street, Wellesley Hills, MA 02481